

BLICKMAN



Tables

Adaptable work and prep surfaces

Includes:

- Instrument Tables
- Utility Tables
- Work Tables




THE BLICKMAN DIFFERENCE

Adjustable Tables

- Our adjustable tables are instrument tables designed to straddle the OR table using either manual or motorized adjustable height mechanisms to accommodate changes in the table position during neuro, cardiac, and orthopedic surgery, as well as varying heights of personnel during general surgery.
- You can customize the tabletop dimensions to fit your differing surface area needs and the leg clearance dimensions to fit your differing OR table dimensions.
- Models are made out of stainless steel and include casters for mobility.
- The motorized table moves along telescoping aluminum lifting columns.

Key Product Features

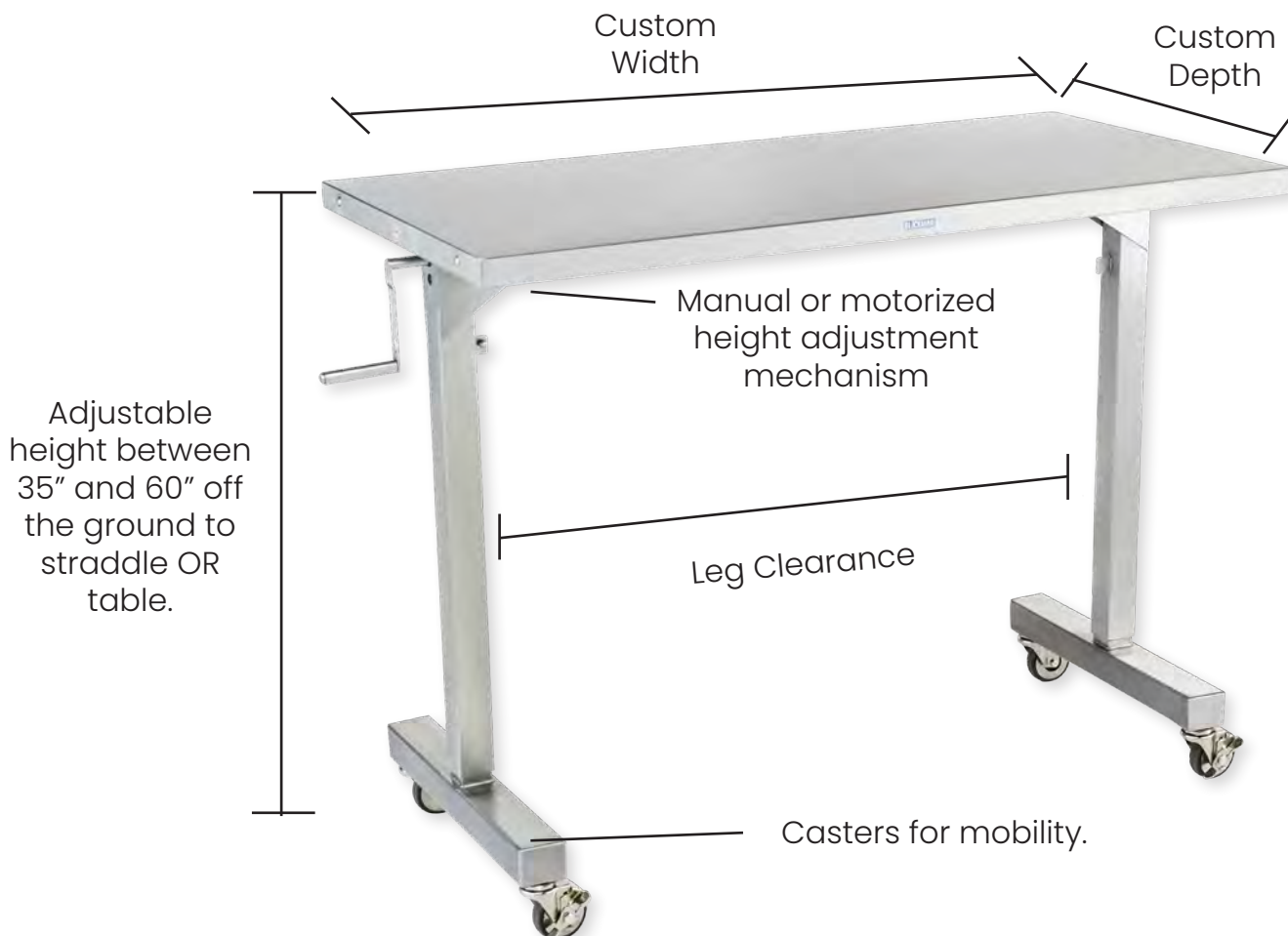

Tables

Adjustable Tables

BACK TABLE

Description: Our specialty back table is designed to provide additional room for extensive surgeries which require a wide array of equipment and instruments. Height adjustment from 46-49" is achievable with either pneumatic assist or manually. The top tier 'hide-away' shelf is angled for easier instrument viewing and holds multiple trays. Tables have sanitary closed tops which help prevent bacterial build up. Sealed reinforcing channels provide rigidity and durability.

Performance Back Table



Manual Back Table



Pneumatic Assist

Model #	Item #	Description	Length	Depth	Height	Overall Adjustable Height	Shipping Weight
BT4830	0364830340	Back Table with 3" Overshelf	48"	30"	34"	46-49"	250 lbs.
BT6030	0366030340	Back Table with 3" Overshelf	60"	30"	34"	46-49"	298 lbs.
BTP6030	036603034P	Back Table Performance with 3" Overshelf	60"	30"	34"	46-49"	360 lbs.
BT7230	0367230340	Back Table with 3" Overshelf	72"	30"	34"	46-49"	345 lbs.
BTP7230	036723034P	Back Table Performance with 3" Overshelf	72"	30"	34"	46-49"	413 lbs.

**All models come standard with 3" casters.*

Custom Dimensions Available Upon Request

Please email us at info@blickman.com


Tables

Adjustable Tables


THE LIBERTY TABLE

Description: Our Liberty table was designed to straddle the OR table with adjustable height mechanisms. It travels effortlessly on antistatic swivel casters with brakes and utilizes a motorized height adjustment mechanism. Use a foot pedal or hand remote to operate the motorized adjustment. Two standard surface areas available.

Customize the surface area dimensions to fit your needs.

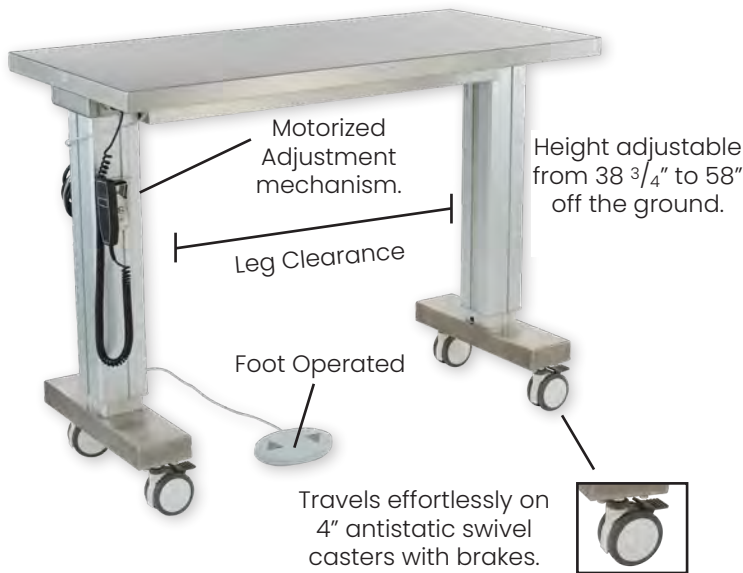

Standard Features

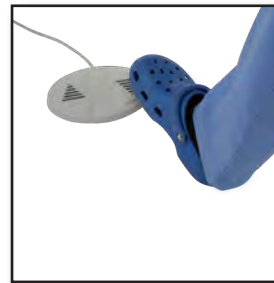

Foot Pedal



Hand Remote

MOTORIZED ADJUSTMENT:

Use the motorized foot operated/hand operated height adjustment mechanism.

Model #	Item #	Description	Width	Depth	Clearance	Shipping Weight
7902SS-M	0157902000	Liberty Table with Adjustable Height 38 3/4-58"	41"	46"	29"	315 lbs.
7901SS-M	0157901000	Liberty Table with Adjustable Height 38 3/4-58"	50 1/2"	24"	38 3/8"	290 lbs.

Custom Dimensions Available Upon Request

Please email us at info@blickman.com

Tables

Adjustable Tables

HYDRAULIC INSTRUMENT TABLE

Description: Our hydraulic instrument table is height adjustable to ensure compatibility with multiple OR table heights and healthcare personnel heights. Adjust its height with a manual hand crank from 30" - 42". It includes one lower shelf and stands stationary on glider feet tips. Select the surface area with dimensions that fit your needs.

Customize the surface area dimensions to fit your needs.



Model #	Item #	Description	Length	Width	Shipping Weight
7835HYD	0177835000	Hydraulic Table with Adjustable Height 30-42"	48"	24"	145 lbs.
7836HYD	0177836030	Hydraulic Table with Adjustable Height 30-42"	60"	30"	188 lbs.

Custom Dimensions Available Upon Request
Please email us at info@blickman.com


Tables

Adjustable Tables

THE GEORGETOWN TABLE

Description: Our Georgetown table was designed to straddle the OR table with a manual hand crank adjustable height mechanism. Height adjustment within ergonomic range between 36" and 56" off the ground. Travels easily on 3" swivel casters with breaks. Two standard surface areas available. Select extended arm design for electrophysiology procedures.

Customize the surface area dimensions to fit your needs.

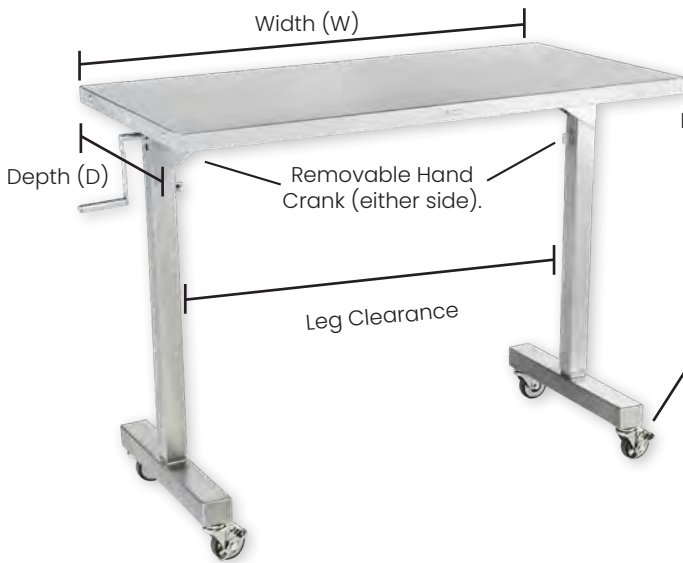

Height adjustable
from 36" to 56" off
the ground.



Travels easily
on four 3" swivel
casters, two with
brakes.

Standard Features




Hand Crank


MANUAL ADJUSTMENT:
Manually adjust the height of this table from
either side with a detachable crank handle.

62" CRT table


This Georgetown table with the
14" extended arm was designed for
optimizing electrophysiology
procedures.

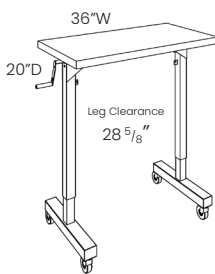


7892SS-0

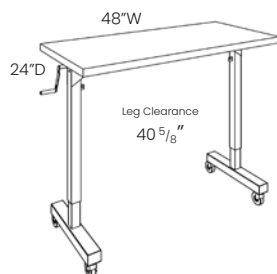

7895SS-0

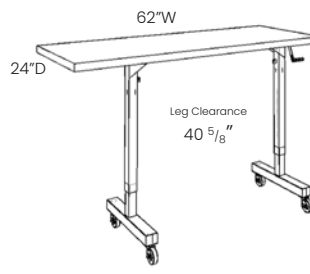

7895SS-0 CRT

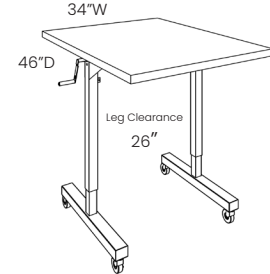

7896SS-0

Name

Georgetown Table

Model #	Item #	Width	Depth	Leg Clearance	Shipping Weight
7892SS-0	0157892000	36"	20"	28 5/8"	135 lbs.
7895SS-0	0157895000	48"	24"	40 5/8"	160 lbs.
7895SS-0 CRT	0157895062	62"	24"	40"	190 lbs.
7896SS-0	0157896000	34"	46"	26"	175 lbs.

Custom Dimensions Available Upon Request

Please email us at info@blickman.com

Tables

Adjustable Tables

OVERBED TABLE

Description: Our overbed table comes in an attractive mark resistant laminate top in a walnut wood grain. It moves up or down with a simple one touch adjustment. Security features include a built-in column lock and four point base. The table travels smoothly on four 2" hooded casters.

Overbed Table Dimensions

Width: 30" x Height: 30" (max) x Depth: 15 1/2"



Model #	Item #	Description	Width	Depth	Height	Shipping Weight
3400	2833400000	Overbed Table	30"	15 1/2"	30" Max	27 lbs.

Construction Specifications

Universal Dimensions	Adjustable heights to fit over OR table or hospital bed
Materials	304 stainless steel Liberty: Telescoping aluminum lifting columns. Overbed Table: Mar resistant laminate (tray), chrome plate steel (column & base)
Gauge	Liberty: 16GA Top; 11GA Sub Top; 14GA Base Georgetown: 18GA Top; 16GA Sub Top, 2" 16GA square tube legs
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish
Casters	Liberty: 4" antistatic swivel casters with brakes Georgetown: 3" swivel casters, two with brakes
Notes	Liberty: Motorized adjustment mechanism Electrical Ratings: 100-120VAC, 15AMP, 60 Hz. Georgetown: Manual adjustment mechanism

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Custom Dimensions Available Upon Request
Please email us at info@blickman.com

THE BLICKMAN DIFFERENCE

Stationary Tables

- We offer a variety of stationary tables to fit any healthcare need. Our stationary tables include instrument tables, utility tables, and general work tables. Some models also come in MR conditional equivalents.
- Our instrument tables are reinforced with H-braces for stability or include an optional lower shelf for storing supplies. We offer standard dimensions or that can be customized to any size. We also have an optional nested set of up to six tables that can be ordered together or separately.
- Our utility tables come in a variety of drawer and shelf configurations to store and protect equipment and can include a workspace so you can store and prep your supplies on one efficient and transportable, setup.
- Our work tables are heavy duty with large surface areas designed for use in sterile and non-sterile preparation.

Key Product Features



Tables

Instrument Tables

NESTED TABLES

Description: Our nested instrument tables are space saving alternatives with a unique design that requires a mere 8 square feet of floor space, but yields 31 square feet of usable work space. All tables travel on 3" casters.



Name	Dimensions (WxHxD)	Model #	Item #	Shipping Weight
Complete Set of 6	8 Sq. Ft. of Floor Space 31 Sq. Ft. of Work Space	782ISS	0137821000	255 lbs.

ORDER SEPARATELY

Largest	48" x 42" x 24"	7822SS	0137822000	56 lbs.
Second	44" x 40" x 22"	7823SS	0137823000	52 lbs.
Third	40" x 38" x 20"	7824SS	0137824000	51 lbs.
Fourth	36" x 36" x 18"	7825SS	0137825000	37 lbs.
Fifth	32" x 34" x 16"	7826SS	0137826000	34 lbs.
Smallest	28" x 32" x 14"	7827SS	0137827000	32 lbs.

MR options available, see MR section for more details



Tables

Instrument Tables

HOWARD & SAWYER TABLES

Description: Our instrument tables offer complete 18 gauge, 304 stainless steel construction. They are lightweight, yet extremely durable. Sanitary closed top help prevents bacterial build up. Sealed reinforcing channels provide rigidity. All tables are 34" tall, which is in the safety zone for reduced back stress, and travel on 3" casters. There is a 16" clearance between shelves. Our 4 most common sizes are listed below, however we offer a wide variety of other sizes that you can order by calling our customer service or visiting our website.

Sawyer Table



Howard Table



Table	Sawyer Table			Howard Table		
Dimensions (WxD)	H-Brace			Shelf		
20" x 16"	7840SS	0127840000	28 lbs.	7830SS	0117830000	32 lbs.
30" x 16"	7839SS	0127839000	33 lbs.	7829SS	0117829000	40 lbs.
33" x 18"	7841SS	0127841000	39 lbs.	7831SS	0117831000	47 lbs.
36" x 20"	7842SS	0127842000	42 lbs.	7832SS	0117832000	52 lbs.
36" x 24"	7843SS	0127843000	48 lbs.	7833SS	0117833000	56 lbs.
48" x 20"	7844SS	0127844000	53 lbs.	7834SS	0117834000	66 lbs.
48" x 24"	7845SS	0127845000	56 lbs.	7835SS	0117835000	69 lbs.
60" x 24"	7846SS	0127846000	178 lbs.	7836SS	0117836000	190 lbs.
72" x 24"	7847SS	0127847000	186 lbs.	7837SS	0117837000	210 lbs.

Options: Custom dimensions, 40"H tables, three sided guardrails, and drawers (see Crescent table design).

MR options available, see MR section for more details



ANGULAR INSTRUMENT TABLE

Description: Choose our angular instrument table to have separate left and right work spaces. Add a lower shelf or drawers for additional storage. This 34" tall table is in the safety zone for reduced back stress, and it travels on 3" casters. Select the surface area with dimensions that fit your needs.



Model #	Item #	Description	Width	Depth	Height	Shipping Weight
7862SSA	0147862000	Angular Instrument Table	60"	15"	34"	120 lbs.
7863SSA	0147863000	Angular Instrument Table	72"	18"	34"	137 lbs.

Custom Dimensions Available Upon Request
Please email us at info@blickman.com

Tables

Utility
Tables

UTILITY TABLES

Description: Our utility tables are popular in the OR as anesthesia tables and throughout healthcare facilities as a general utility table. Our superior construction adds a supporting channel underneath the top shelf to enhance its strength. All of our utility tables are 34" tall, which is in the safety zone for reduced back stress, and travel easily on 3" casters.

Shelf Configuration Includes:

- ✓ Two shelves, 17" and 34" off the ground
- ✓ Four sided guardrails for protection

Drawer Configurations Include:

- ✓ A positive stop to prevent accidental drawer removal
- ✓ Ball bearing slides for easy sliding drawers
- ✓ Three sided guardrails for protection (optional four sided guardrails are available)

All of our utility tables are the same dimensions:

Width: 20" x Height: 34" x Depth: 16"



Lower Shelf Clearance:
16"

Lower Shelf Clearance:
11"

Lower Shelf Clearance:
5 1/2"

Lower Shelf Clearance:
4 1/4"

Lowest Drawer Clearance:
11 1/2"

# of Drawers	0	1	2	3	4
Model #	7852SS	7850SS	7754SS	7857SS	7859SS
Item #	0257852000	0217850000	0227754000	0237857000	0247859000
Shipping Weight	32 lbs.	46 lbs.	59 lbs.	72 lbs.	82 lbs.

Options: Feet: glides, rubber tips, or larger diameter casters, four sided guard rails on drawer models

Tables

CRESCENT TABLES

Description: Our crescent table has the same construction as the Howard and Sawyer Instrument tables with the addition of drawers and guard rails. This table includes a large work surface area (36"x 20"), two 4" deep drawers and one optional lower shelf 16" off the ground. This model is 34" tall, which is in the safety zone for reduced back stress, and travels easily on 3" casters.

Customize the surface area dimensions of these models to fit your needs.



Model #	Item #	Description	Width	Depth	Height	Shipping Weight
7856SS	0317856000	Crescent Table (H-Brace)	36"	20"	34"	78 lbs.
7856SS-SH	0317856003	Crescent Table (Shelf)	36"	20"	34"	81 lbs.

Options: Feet: glides, rubber tips, or larger diameter casters

WORK TABLES

Description: Our work tables are designed for use as prep tables. They are heavy duty with large surface areas designed for use in sterile and non-sterile environments. The smallest table dimensions use four legs, and all other larger dimensions use six legs. All of our work tables are 36" tall, which is in the safety zone for reduced back stress, and come standard with bullet feet. Rubber tips, casters, and gliders are available as alternative options.

Customize the surface area dimensions of these models to fit your needs.



Model #	Item #	Description	Width	Depth	Height	Shipping Weight
H-Brace						
CSDT6830	0336830000	Work Table	68"	30"	36"	210 lbs.
CSDT8036	0338036000	Work Table	80"	36"	36"	300 lbs.
CSDT9244	0339244000	Work Table	92"	44"	36"	390 lbs.
Shelf						
CSDT6830S1	0336830100	Work Table	68"	30"	36"	260 lbs.
CSDT8036S1	0338036100	Work Table	80"	36"	36"	355 lbs.
CSDT9244S1	0339244100	Work Table	92"	44"	36"	435 lbs.

Options: Feet: glides, rubber tips, or larger diameter casters

Custom Designs Available Upon Request
Please email us at info@blickman.com

STATIONARY TABLE OPTIONS

Drawers

Description: For added storage, add drawers beneath any tabletop.



**contact us for available drawers and sizing*

Guardrails

Description: For added placement security, add guardrails to any tabletop.



Custom Designs Available Upon Request
Please email us at info@blickman.com

Tables

Options

STATIONARY TABLE OPTIONS

Feet Options

Rubber Tips | Casters | Gliders



Custom Designs Available Upon Request

Please email us at info@blickman.com

Tables

STATIONARY TABLE
Construction
Specifications

Construction Specification

Universal Dimensions	Custom dimensions available on all tables. Please email us at info@blickman.com .
Materials	304 stainless steel
Gauge	Instrument Tables: 1 1/4" tube legs, 18GA table top Utility Tables: 1" Diameter 18GA tube legs, 18GA top & body Work Tables: 1 5/8", 16GA tube legs, 16GA table top & sub top
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish
Casters	3" Casters
Notes	MR conditional tables available where indicated

Custom Dimensions/Designs Available Upon Request

Please email us at info@blickman.com

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Toll Free: (800) 247-5070
Local: (973) 330-0557

BLICKMAN

www.blickman.com
info@blickman.com

1099 Wall Street West, Suite 250A, Lyndhurst, NJ 07071