

# BLICKMAN



## ***MR Conditional Room***

*Broadest line of MR equipment*

**Includes:**

- Exam Stools
- Hampers
- Foot Stools
- Stands
- Tables



# THE BLICKMAN DIFFERENCE

## ***How we make our products “MR Conditional”***

We manufacture products primarily from type 304 stainless steel, which is regarded as non-magnetic in its original, annealed condition. Forming the metal tends to bring out the magnetizing properties, so we demagnetize the metal via solution annealing and a Nitrogen quenching process. We also make certain to employ, where appropriate, non-metal parts and components such as fasteners, hinges, and casters.

## ***Broadest Line of MR Equipment in the Healthcare Market***

- Padded stools on rubber tips or on casters
- Single platform or double platform foot stools with or without handrails
- Rolling hampers, with or without plastic lids
- Rolling IV stands, mayo stands, solution stands, and kickbuckets
- Instrument tables and utility tables in a variety of sizes with or without shelves/casters

## ***Fully compliant with the New ASTM Standards***

- All Blickman products are labeled MR conditional or *NOT* MR conditional.
- All Blickman MR Equipment is suitable for use throughout the entire MR suite, including anywhere within the vicinity of 3.0 Tesla scanners.



## ***Backed by our Limited Lifetime Warranty***

Since the late 1990s, our products, legendary for high quality and supported by an industry-unique Limited Lifetime Warranty, have performed without incident. Our exceptional record of safety in the imaging industry is a direct result of our consistent adherence to guidelines set by the Center for Devices & Radiological Health (CDRH) - Magnetic Resonance Group, part of the U.S. Food and Drug Administration (FDA).



***For details, please visit:***

***<https://www.blickman.com/limited-lifetime-warranty>***

# MR Conditional

## Exam Stools

### EXAM STOOL

**Description:** Our MR conditional exam stool is industry standard, fully welded, stainless steel construction. Sit comfortably on 15" diameter, 3" thick padded seats. Choose between stationary feet tips or casters for mobility. Height is 21" on feet tips and 23" on MR casters.



Name	MR Conditional Stool	MR Conditional Stool with Casters
Item #	7445MR	7445MR-C
Model #	1027445000	1027445001
Shipping Weight	14 lbs.	16 lbs.

# MR Conditional Hamper

## HAMPER

**Description:** Our MR conditional round stainless steel hamper travels easily on MR casters and is compatible with our 18" nylon hamper bags and our water soluble hamper bags.



---

Name	MR Conditional Hamper
Model #	7772MR
Item #	0927772400
Shipping Weight	20 lbs.


---

# MR Conditional

## Foot Stools

### KENT FOOT STOOLS

Our MR conditional Kent foot stools include one step 8" off the ground. They are industry standard, fully welded, stainless steel construction with non-slip rubber treads and feet tips. 18" and 30" width sizes are available. Include a 41" tall handrail for additional support.

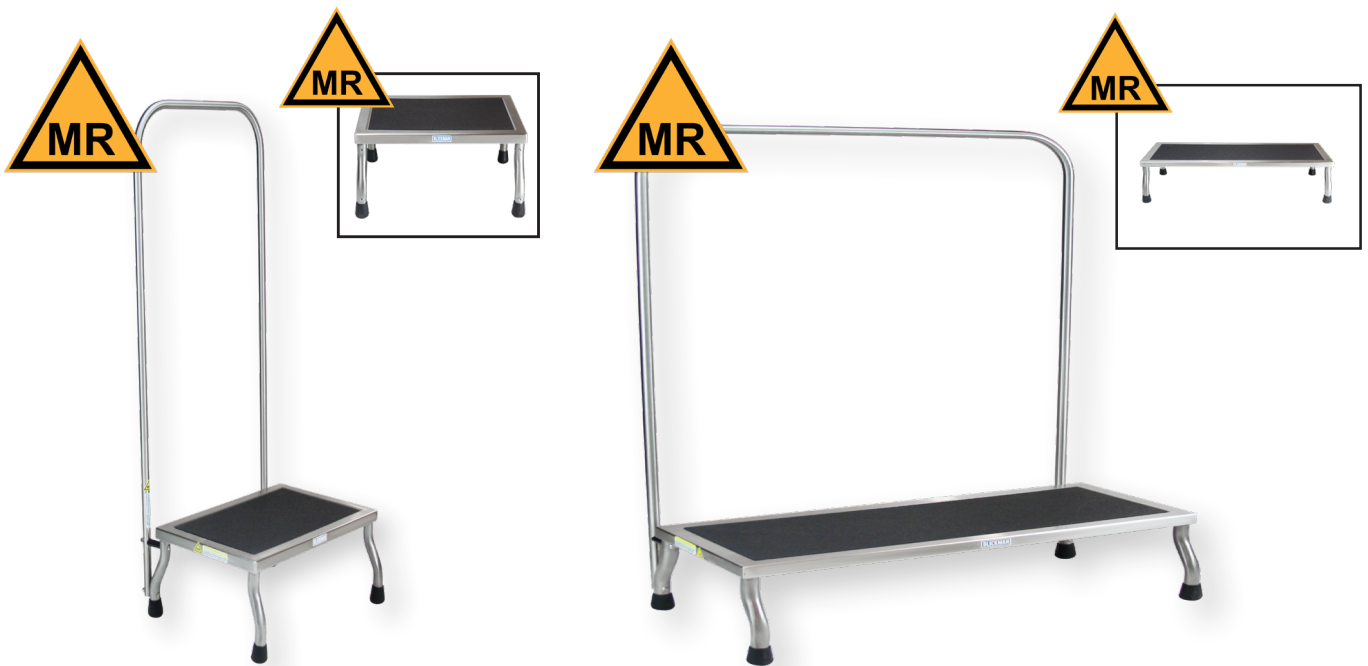

Name	MR Conditional 18" Kent Stool	MR Conditional 30" Kent Stool
	<b>With Handrail (41" high)</b>	
Dimensions (WxHxD)	18" x 41" x 12"	30" x 41" x 12"
Model #	7757MR-HR	7761MR-HR
Item #	1017757100	1017761100
Shipping Weight	20 lbs.*	30 lbs.*
	<b>Without Handrail</b>	
Dimensions (WxHxD)	18" x 8" x 12"	30" x 8" x 12"
Model #	7757MR	7761MR
Item #	1017757000	1017761000
Shipping Weight	11 lbs.	20 lbs.

# MR Conditional

Two Step

## TWO STEP FOOT STOOLS

**Description:** Our MR conditional two step foot stools are industry standard, fully welded, stainless steel construction with H-braces for additional reinforcement, non-slip rubber treads for secure footing, and feet tips for firm floor placement.



**Step Surface Areas**  
 Top Step: 36" x 12"  
 Bottom Step: 36" x 24"

**Heights of Steps**  
 Top Step: 15 3/4"  
 Bottom Step: 8"

Name	Two Step Foot Stool
	<b>With Handrails (41" high)</b>
Model #	7764MR-HR
Item #	1012041096
Shipping Weight	65 lbs.

### Two Step Foot Stool with Side Handrails



**Step Surface Areas**  
 Top Step: 24" x 8"  
 Bottom Step: 24" x 16"

**Heights of Steps**  
 Top Step: 15 3/4"  
 Bottom Step: 8"

Name	The Donnelly
	<b>With Handrails (41" high)</b>
Model #	7763MR-HR
Item #	1017763100
Shipping Weight	40 lbs.*
	<b>Without Handrail</b>
Model #	7763MR
Item #	1017763000
Shipping Weight	27 lbs.

### The Donnelly with Handrail

\*Shipping weight with an asterisk is the combined shipping weight of separate parts

# MR Conditional Stands

## NEWARK MAYO STAND

**Description:** Our stainless steel MR conditional Mayo stand has a 12 5/8" x 19 1/8" stainless steel tray and travels easily on MR casters. It is height adjustable from 39 1/2" to 62" using the Tru-Loc™ friction knob.



Name	MR Conditional Newark Mayo Stand
Model #	8841MR
Item #	0638841200
Shipping Weight	36 lbs.

## IV STAND

**Description:** Our 5 leg MR conditional IV stand consists of a stainless steel pole on a 24" diameter, polished diecast aluminum base with Wall Saver™ caps and MR casters. It is height adjustable from 74" to 110" via a twist lock adjustment mechanism. Choose from either 2- or 4-hooks or contact customer service for 6 or 8 (rake) hooks.



Name	MR Conditional IV Stands
<b># of Hooks</b>	<b>2 Hooks</b>
Model #	7792MR-2
Item #	0537792002
Shipping Weight	16 lbs.
<b># of Hooks</b>	<b>4 Hooks</b>
Model #	7792MR-4
Item #	0537792004
Shipping Weight	17 lbs.

### MR Push Handle for MR IV Stands



Name	Model #	Item #	Ship Weight
Push Handle	1381MR	3537426350	5 lbs.

# MR Conditional

Stands

## SOLUTION STANDS

**Description:** We offer two types of stainless steel MR conditional solution stands: a single basin model and a double basin model. Our solution stands are 33" tall offering compatibility with surgical table heights and reduced back stress. Choose between a lower shelf or an H-brace for added storage and support.



### Single Basin

**Dimensions:**

Width: 15" x Height: 33" x Depth: 15"



### Double Basin

**Dimensions:**

Width: 29" x Height: 33" x Depth: 15"

Name	MR Conditional Single Basin SS	MR Conditional Double Basin SS
<b>Supporting Frame</b>		<b>Shelf</b>
Model #	7807MR	7808MR
Item #	0717807066	0727808066
Shipping Weight	20 lbs.	33 lbs.
<b>Supporting Frame</b>		<b>H-Brace</b>
Model #	7807MR-HB	7808MR-HB
Item #	0717807166	0727808166
Shipping Weight	19 lbs.	24 lbs.

Note: SS denotes solution stand

## KICKBUCKET

**Description:** Our stainless steel MR conditional kickbucket consists of a continuous encircling ubber bumper, which helps prevent tipping, wall scuff marks, and noise. A stainless steel cross brace provides excellent protection and a secure hold of a heavy gauge 13 qt. pail, and MR casters for mobility.



Name	Model #	Item #	Ship Weight
MR Conditional Lenox Kickbucket with Pail	7766MR	0817766066	15 lbs.



# MR Conditional

## Tables

### NESTED TABLES

**Description:** Our MR conditional nested instrument tables are space saving alternatives with a unique design that requires a mere 8 square feet of floor space, but yields 31 square feet of usable work space. All nested tables travel on MR casters.



Name	Dimensions (WxHxD)	Model #	Item #	Ship Weight
Complete Set of 6	8 Sq. Ft. of Floor Space 31 Sq. Ft. of Work Space	7821MR	0197821000	255 lbs.

**ORDER SEPARATELY**

Largest	48" x 42" x 24"	7822MR	0197822100	56 lbs.
Second	44" x 40" x 22"	7823MR	0197823100	52 lbs.
Third	40" x 38" x 20"	7824MR	0197824100	51 lbs.
Fourth	36" x 36" x 18"	7825MR	0197825100	37 lbs.
Fifth	32" x 34" x 16"	7826MR	0197826100	34 lbs.
Smallest	28" x 32" x 14"	7827MR	0197827100	32 lbs.

# MR Conditional

Tables

## INSTRUMENT TABLES

**Description:** Our MR conditional instrument tables offer complete 18 gauge, 304 stainless steel construction. A sanitary closed top helps prevent bacterial build up. Sealed reinforcing channels provide increased rigidity. All tables are at an ergonomic height of 34" tall with 1-1/4" diameter legs that are fully welded for strength and stability and travel easily on MR casters. There is a 16" clearance between shelves. Our 4 most common sizes are listed below, however you can order a customized size by calling customer service or visiting our website.



Surface Area	MR Conditional Sawyer Table	MR Conditional Howard Table
<b>20" x 16"</b>		
Model #	7840MR	7830MR
Item #	0197840100	0187830100
Shipping Weight	28 lbs.	32 lbs.
<b>36" x 24"</b>		
Model #	7843MR	7833MR
Item #	0197843100	0187833100
Shipping Weight	48 lbs.	56 lbs.
<b>48" x 24"</b>		
Model #	7845MR	7835MR
Item #	0197845100	0187835100
Shipping Weight	56 lbs.	69 lbs.
<b>60" x 24"</b>		
Model #	7846MR	7836MR
Item #	0197846100	0187836100
Shipping Weight	189 lbs.	200 lbs.

**Custom Dimensions Available Upon Request**

Please email us at [info@blickman.com](mailto:info@blickman.com)

# MR Conditional

## Exam Stools

### UTILITY CARTS

**Description:** Our MR conditional utility carts are made out of industry standard, fully welded, stainless steel construction. They feature two to three easy-access storage shelves and two ergonomic push handles. Optional IV pole attachments available for every model.



Name	Medium Duty End Rails Utility Cart	Medium Duty Side Rails Utility Cart	Heavy Duty Utility Cart
MR Conditional Standard Models			
Item #	7536SS-MR	7535SS-MR	7537SS-MR
Model #	2497536000	2497535000	2497537000
Shipping Weight	40 lbs.	47 lbs.	43 lbs.
MR Conditional Models with IV Attachments			
Item #	7536SS-IV-MR	7535SS-IV-MR	7537SS-IV-MR
Model #	2497536001	2497535001	2497537001
Shipping Weight	45 lbs.	52 lbs.	48 lbs.

# MR Conditional

## Specifications

### Construction Specifications

All products in this Blickman MR Conditional Room brochure meet the revised “MR Conditional” definition as stated in ASTM standard F2503.

The configuration and design safety of the components of these Blickman products in the MR environment were verified in accordance with ASTM standard F2052-06: “Measurement of Magnetically-induced Displacement Force in the MR Environment.”

The “paramagnetic” material property constraints of these Blickman products were determined through testing conducted by an independent lab in accordance with ASTM standard A773/A773M-01, employing a computer automated, soft magnetic, hysteresis graph system.

A low- $\mu$  permeability indicator gage, in compliance with ASTM standard A342/A342M, is used for final factory acceptance of these Blickman products.

Identification of these Blickman products conforms to the ASTM standard F2503-05: “Marking Medical Devices and Other Items for Safety in the MR Environment.”

### Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at [info@blickman.com](mailto:info@blickman.com).

Toll Free: (800) 247-5070  
Local: (973) 330-0557

**BLICKMAN**

[www.blickman.com](http://www.blickman.com)  
[info@blickman.com](mailto:info@blickman.com)

1099 Wall Street West, Suite 250A, Lyndhurst, NJ 07071