

BLICKMAN



Mayo Stands

Maximum utility in minimal space

Includes:



- Chrome Mayo stands
- Stainless steel Mayo stands



THE BLICKMAN DIFFERENCE

- Our Mayo stands have evolved to meet the needs of the marketplace for maximum utility in an increasingly crowded surgical environment.
- We have an array of options with varying material, bases, shapes, and adjustment mechanisms to suit your needs.
- All of our Mayo stands are mobile and use 2" roller blade, hard rubber, or twin wheel plastic hooded casters.

Key Product Features

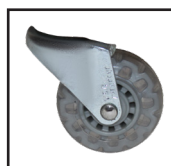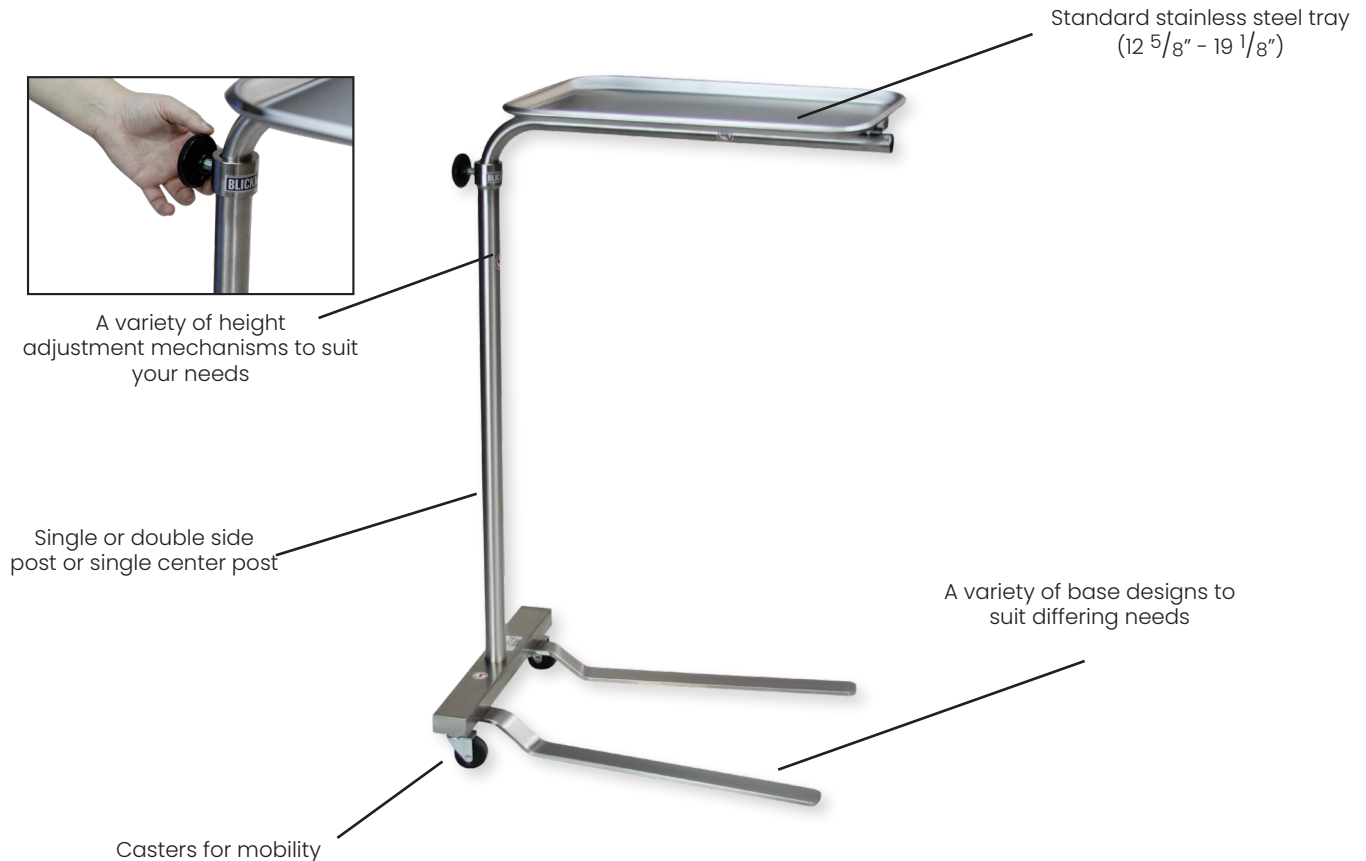

2" Rollerblade Casters



2" Hard Rubber Casters



2" Twin Wheel Hooded Casters

Mayo Stands

Chrome

CHROME MAYO STANDS

Description: Each of our chrome Mayo stands contain a 12 5/8" x 19 1/8" stainless steel tray and are height adjustable using the Tru-Loc™ friction knob. We offer different base and pole set ups to meet your varying needs.

Middle Pole

Base: Five leg cast aluminum base with wall savers™ and 2" twin wheel hooded casters



Double Pole

Base: Chrome loop base with two 2" twin wheel hooded back casters



Single Pole

Base: Two pronged chrome base with single hard rubber back casters



| Name | Middle Pole | Double Pole | Single Pole |
|-------------------|---------------|---------------|---------------|
| Height Adjustment | 29 1/2" - 49" | 34 1/2" - 54" | 34 1/2" - 54" |
| Model # | 1501 | 1510 | 1515 |
| Item # | 0661501000 | 0661510000 | 0661515000 |
| Shipping Weight | 14 lbs. | 18 lbs. | 12 lbs. |

Mayo Stands Stainless Steel

NEWARK MAYO STANDS

Description: Each of our stainless steel Mayo stands contain a 12 5/8" x 19 1/8" stainless steel tray and are height adjustable using the Tru-Loc™ friction knob. We offer different base set ups to meet your varying needs.

Lightweight

Base: Two pronged stainless steel base with two hard rubber back casters



Standard



Base: Two pronged stainless steel base with two roller blade back casters



Four Caster

Base: Two pronged stainless steel base with four roller blade casters



| Name | Lightweight |  Standard |  Four Caster |
|-------------------|---------------|--|---|
| Height Adjustment | 37 1/4" - 64" | 39 1/2" - 62" | 39 1/2" - 62" |
| Model # | 1520SS | 8841SS | 8841MR |
| Item # | 0671520000 | 0638841000 | 0638841200 |
| Shipping Weight | 26 lbs. | 36 lbs. | 36 lbs. |




Mayo Stands Stainless Steel

BENJAMIN MAYO STAND

Description: Our Benjamin Mayo stand design is made for ease and durability. Our state of the art double post design minimizes lateral movement and increases overall stability. We offer this design in three different tray sizes.



Via a foot operated mechanism, you can adjust the height from 39 1/2" to 62" high off the ground. A guard over the foot prevents accidental release.

| Name | Benjamin Mayo Stands | |
|------------------|-------------------------------------|---|
| Tray Size | 12 5/8" x 19 1/8" (Standard) |  |
| Model # | 8867SS | |
| Item # | 0668867000 | |
| Shipping Weight | 40 lbs. | |
| Tray Size | 16" x 21" |  |
| Model # | 8869SS | |
| Item # | 0668869000 | |
| Shipping Weight | 41 lbs. | |
| Tray Size | 20" x 25" |  |
| Model # | 8871SS | |
| Item # | 0668871000 | |
| Shipping Weight | 43 lbs. | |

Mayo Stands

Construction Specifications

Construction Specifications

| | |
|----------------------|---|
| Universal Dimensions | Standard Tray Size: Width: 12 5/8" x Depth: 19 1/8" |
| Materials | 304 stainless steel Chrome |
| Gauge | All contain stainless steel trays: 20GA Benjamin Mayo Stands: 3/4" diameter, 18GA tubes, 1/8" base Newark Mayo Stands: 1520SS: 1 1/4" diameter, 18GA tube, weighted base 8841SS/8848SS: 1" diameter, 16GA tube, 3/8" base |
| Welding | All welding by heliarc method and exposed welds are ground smooth |
| Finish | Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish |
| Casters | 2" plastic hooded twin wheel, hard rubber, or roller blade casters |
| Notes | <ul style="list-style-type: none">MR conditional Mayo stands available where indicated |

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Toll Free: (800) 247-5070
Local: (973) 330-0557

BLICKMAN

www.blickman.com
info@blickman.com

1099 Wall Street West, Suite 250A, Lyndhurst, NJ 07071