

# BLICKMAN



## ***Scrub Sinks***

*Hands free multi-station operation*

**Includes:**

- Windsor Scrub Sinks
- Lodi Scrub Sinks
- Accessories




# THE BLICKMAN DIFFERENCE

## Scrub Sinks

- Our two scrub sink models, the Windsor and the Lodi, are each available with single, dual, and triple stations. A four-station unit is also available in the Windsor model.
- They have deep sink bowls with sound-deadened enclosed double walled panels that are angled in to divert water splashes along with spouts positioned toward the rear of the sink. Standard gooseneck faucets included.
- All operation options are hands free. The Lodi scrub sinks have infrared water control and knee-action water and soap dispensers, and the Windsor scrub sinks come standard with knee-action control, but they are also available in foot, elbow, and infrared sensor controls.

## Key Product Features

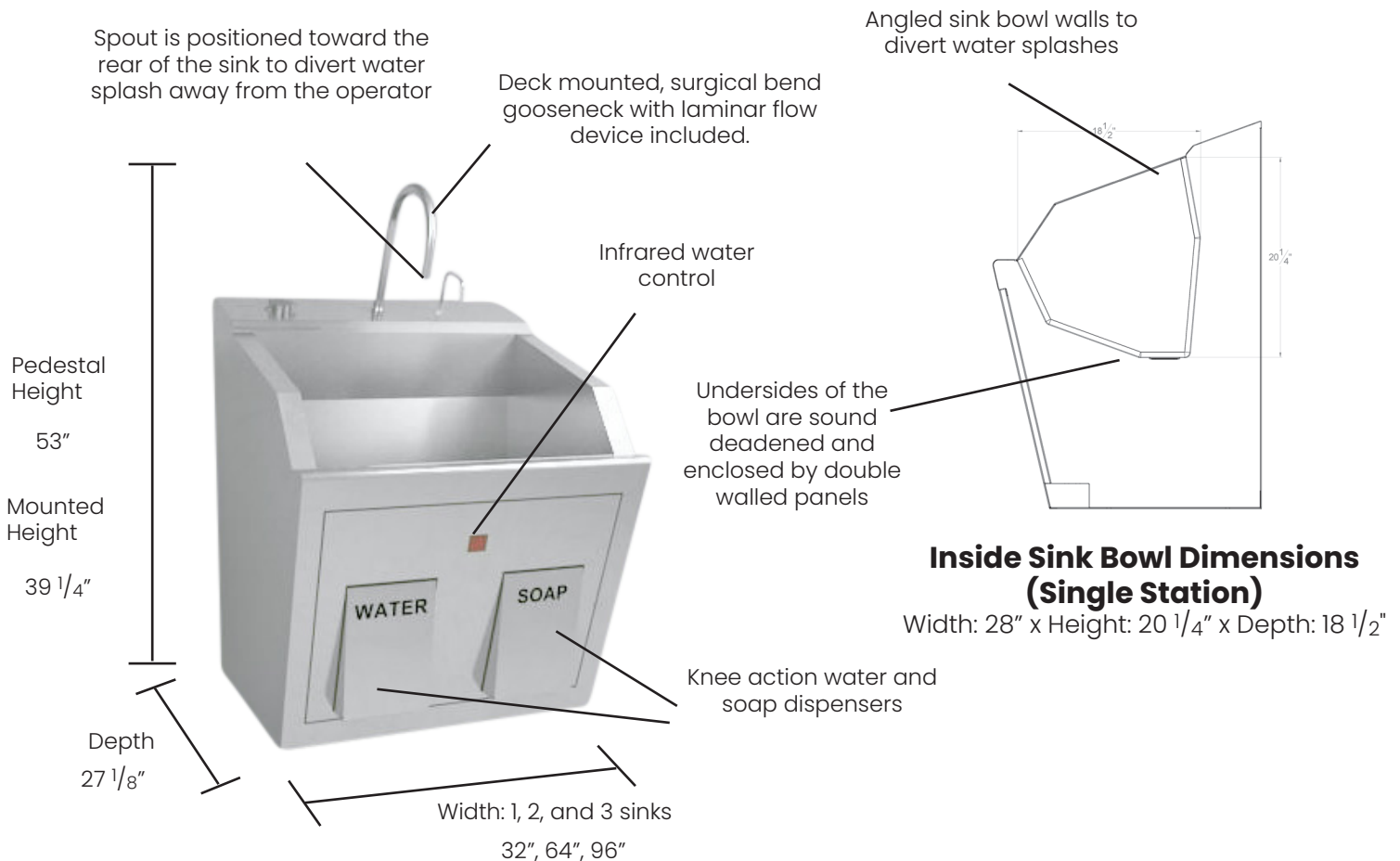

# Scrub Sinks

Windsor

## WINDSOR SCRUB SINKS

**Description:** Our Windsor scrub sinks are all wall mounted and utilize one undivided sink bowl. Choose from our four different hands free operations: knee, foot, elbow, or infrared. Standard gooseneck faucets included.

### Inside Sink Bowl Dimensions

Width: 21 3/4" x Height: 8" x Depth: 14 1/2"



**Width**  
25"



**Width**  
47"



**Width**  
69"



**Width**  
91"

Name	1 Station	2 Stations	3 Stations	4 Stations
<b>Knee Action</b>				
Model #	7878SS-KA	7879SS-KA	7880SS-KA	7881SS-KA
Item #	1317878000	1317879000	1317880000	1317881000
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.
<b>Foot Action</b>				
Model #	7878SS-FO	7879SS-FO	7880SS-FO	7881SS-FO
Item #	1317878003	1317879003	1317880003	1317881003
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.
<b>Elbow Action</b>				
Model #	7878SS-EA	7879SS-EA	7880SS-EA	7881SS-EA
Item #	1317878004	1317879004	1317880004	1317881004
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.
<b>Infrared</b>				
Model #	7878SS-IF	7879SS-IF	7880SS-IF	7881SS-IF
Item #	1317878001	1317879001	1317880001	1317881001
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.

# Scrub Sinks

Lodi

## LODI SCRUB SINKS

**Description:** Our Lodi scrub sinks come standard with both infrared water control and knee-action water and soap dispensers. All units are available in either pedestal or wall-mounted configurations. Standard gooseneck faucets included.



**Width**  
32"



**Width**  
64"



**Width**  
96"



**Height**  
53"



**Height**  
39 1/4"

Name	1 Station	2 Stations	3 Stations
<b>Pedestal</b>			
Model #	9881SSP	9882SSP	9883SSP
Item #	1339881P00	1339882P00	1339883P00
Shipping Weight	275 lbs.	420 lbs.	600 lbs.
<b>Wall Mounted</b>			
Model #	9881SSW	9882SSW	9883SSW
Item #	1339881W00	1339882W00	1339883W00
Shipping Weight	215 lbs.	360 lbs.	570 lbs.

# Scrub Sinks

Accessories

## LODI SCRUB SINK ACCESSORIES

- Factory Installed Options -



**Digital Timer**



**Eyewash**

- For Wall Mount Installation -



**Chair Carrier**  
For Wall Mount Installation

# of Stations	One	Two	Three
Part #	4301224751	4301224752	4301224753
Shipping Weight	47 lbs.	119 lbs.	141 lbs.

# Scrub Sinks

## Construction Specifications

### Construction Specifications

Standard Dimensions **Lodi:** (W) 32", 64", 96" x (D) 27 1/4" x (H): 39 1/4", 53"  
Internal Sink (Single Station): 28" x 18 1/2" x 20 1/4"  
**Windsor:** (W) 25", 47", 69", 91" x (D) 18 1/2" x (H): 15 3/4"  
Internal Sink: 21 3/4" x 8" x 14 1/2"

Materials/  
Gauge 304 stainless steel  
Lodi: Mix of 16GA and 14GA  
Windsor: 16 GA  
Double/Triple Station Lodi Scrub Sinks: Plastic Dividers

Welding All welding by heliarc method and exposed welds are ground smooth

Finish Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish

Electrical Available in standard 100/120 VAC input or 220/240 VAC  
Units are FDA registered and conform to standards for Safety Medical Electrical Equipment UL 606011 (US) and CSA C22.2 #601.1M90 (Canada)  
Units are cord connected with a hospital grade cap

Casters N/A

### Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at [info@blickman.com](mailto:info@blickman.com).

Toll Free: (800) 247-5070  
Local: (973) 330-0557

**BLICKMAN**

[www.blickman.com](http://www.blickman.com)  
[info@blickman.com](mailto:info@blickman.com)

1099 Wall Street West, Suite 250A, Lyndhurst, NJ 07071