

BLICKMAN



Utility Carts

Convenient transport for medical equipment

Includes:


- Light Duty
- Medium Duty
- Heavy Duty



THE BLICKMAN DIFFERENCE

- Our utility carts are ideal for transporting medical equipment which needs to be used in various areas in your medical facility. We offer light duty, medium duty, and heavy duty options to suit all of your needs.
- Common features include a push handle, open construction for easy access, guardrails for securing contents, and 4" swivel casters for easy maneuverability.
- They contain two or three solid stainless steel shelves to provide ample storage space.
- The stainless steel construction allows for easy cleaning and sanitizing, while resisting corrosion and stains.
- We also offer a folding utility cart with wire baskets for any application such as movement of mail, patient records, and laboratory materials.

Key Product Features

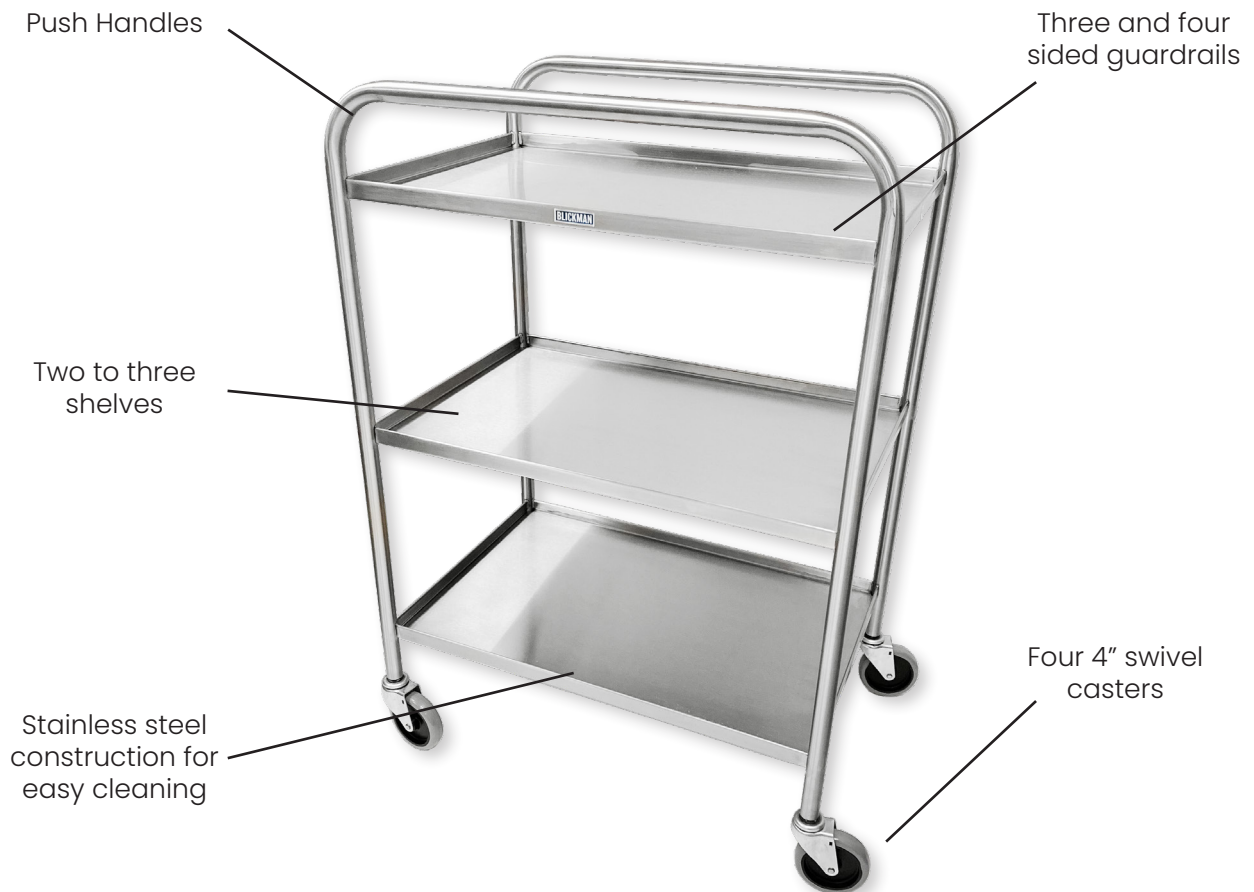

Utility Carts

Light Duty

FOLDING UTILITY CART

Description: Our folding utility cart is versatile enough for any application such as movement of mail, patient records, and laboratory materials. Choose any combination of 6" or 12" deep wire baskets to complete your cart. The metal frame is designed for easy assembly and folds for storage. Rubber stoppers support and protect the wire racks and provide a quiet ride on four, 3" hard rubber casters.

Folding Utility Cart Dimensions:

Width: 29 1/4" x Height: 38 1/2" x Depth: 19 1/2"



6" Wire Basket (2441)

Item #: 2422441000



12" Wire Basket (2442)

Item #: 2422442000

Name	Model #	Item #	Shipping Weight
Kit with 6" + 12" Baskets	2440 KIT	24224400KT	33 lbs.
Kit with Two 6" Baskets	2440 6 KIT	24224401KT	31 lbs.
Kit with Two 12" Baskets	2440 12 KIT	24224402KT	35 lbs.

ORDER SEPARATELY

Folding Frame	2440	2422440000	21 lbs.
6" Wire Basket	2441	2422441000	5 lbs.
12" Wire Basket	2442	2422442000	6 lbs.

Utility Carts

Light Duty

LIGHT DUTY CARTS

Description: Our versatile light duty cart easily moves records, supplies and equipment throughout your facility. It's made of 20GA stainless steel, has three solid shelves, one push handle with protective rubber bumpers on each corner, travels on four 4" swivel casters and has a weight capacity of 300 lbs. evenly distributed. The top and middle shelves have three edges with 1-1/4" sides and a front open edge to slide equipment in and out with ease. The bottom shelf has four open edges for greater accessibility to equipment stored in a lower position. You can also choose to add additional guardrails or an IV pole attachment.

Our light duty carts have the same dimensions:

Width: 28 1/4" x Height: 35 3/4" x Depth: 16"



Top Shelf: Three Sided Shelf with Open Front Edge
Middle Shelf: Three Sided Shelf with Open Front Edge
Bottom Shelf: Open Shelf



Top Shelf: Three Sided Shelf with Open Front Edge
Middle Shelf: Three Sided Shelf with Front Guardrail
Bottom Shelf: Shelf with Four Guardrails



Top Shelf: Three Sided Shelf with Front Guardrail
Middle Shelf: Three Sided Shelf with Front Guardrail
Bottom Shelf: Shelf with Four Guardrails

Name	Light Duty Utility Cart	Light Duty Utility Cart with 5 Guardrails	Light Duty Utility Cart with 6 Guardrails
Model #	7534SS	7534SS-GR	7534SS-GR4
Item #	2427534002	2427534000	2427534004
Shipping Weight	30 lbs.	30 lbs.	30 lbs.

Options: IV pole attachments.

Custom Designs Available Upon Request
 Please email us at info@blickman.com

Utility Carts

Medium/Heavy Duty

MEDIUM DUTY CARTS

Description: Our versatile medium duty carts easily move records, supplies and equipment throughout your facility. They are made of 16GA stainless steel, have three solid shelves with lipped edges on all four sides, travel on four 4" swivel casters and have a weight capacity of 500 lbs. evenly distributed. Choose between our design that includes four sided shelves and width-wise handles or our design that includes open shelving and depth-wise handles. Optional IV pole attachment is available.



Dimensions: 27 1/2" x 40" x 18 1/8"



Dimensions: 25" x 35 3/4" x 16 1/2"

Name	Side Rails Medium Duty Cart	End Rails Medium Duty Cart
Model #	7535SS	7536SS
Item #	2427535000	2427536000
Shipping Weight	47 lbs.	40 lbs.

HEAVY DUTY CART

Description: Our heavy duty cart with angled push handles is made of 16GA stainless steel, has two solid 4-sided shelves, travels on four 5" swivel casters and has a weight capacity of 700 lbs. evenly distributed. Optional IV pole attachment is available.



Dimensions: 41 3/4" x 40" x 19 1/4"

Name	Heavy Duty Cart
Model #	7537SS
Item #	2427537000
Shipping Weight	43 lbs.

Custom Designs Available Upon Request

Please email us at info@blickman.com

MR options available, see MR section for more details.



Utility Carts

Accessories

IV POLE ATTACHMENTS

Description: Our IV pole attachment can fit any Blickman utility cart with ease.



Name	7534SS with IV pole attachment	7535SS with IV pole attachment	7536SS with IV pole attachment	7537SS with IV pole attachment
Model #	7534SS-IV	7535SS-IV	7536SS-IV	7537SS-IV
Item #	2427534001	2427535001	2427536001	2427537001
Shipping Weight	35 lbs.	52 lbs.	45 lbs.	48 lbs.

GUARDRAILS & SIDING

Contact Blickman Customer Service about adding additional siding or guardrails.



Open Shelf



Lipped Shelf



Shelf with 4 Guardrails



3 Sided Shelf

Custom Designs Available Upon Request

Please email us at info@blickman.com

MR options available, see MR section for more details.



Utility Carts

Construction Specifications

Construction Specifications

Dimensions	Light Duty: 28 1/4" W x 35 3/4" H x 16" D Medium Duty: 27 1/2" W x 40" H x 18 1/8" D (Side Rails) 25" W x 35 3/4" H x 16 1/2" D (End Rails) Heavy Duty: 41 3/4" W x 40" H x 19 1/4" D
Materials	Stainless steel Rubber bumper handles on 7534SS Carts
Gauge	Light Duty: 20GA & 18GA Medium/Heavy Duty: 16GA
Welding	All welding by spot and heliarc methods and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with a #4 polish to a satin finish
Casters	Light/Medium Duty: 4" Swivel Casters Heavy Duty: 5" Swivel Casters

Custom Designs Available Upon Request
Please email us at info@blickman.com

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Toll Free: (800) 247-5070
Local: (973) 330-0557

BLICKMAN

www.blickman.com
info@blickman.com

1099 Wall Street West, Suite 250A, Lyndhurst, NJ 07071