

Optimal infection control

16



Includes:

- Decontamination Stations
- Faucet Options
- Accessories

THE BLICKMAN DIFFERENCE

Decontamination Stations

- Our decontamination stations provide optimal infection control in today's high risk healthcare environments. Our durable stainless steel decontamination stations meet the stringent Association for the Advancement of Medical Instrumentation (AAMI) guidelines for maximum sterility and utility.
- Our decontamination stations have an ergonomic working height for reduced back stress, symmetrical design for either left or right work flow, 18" counter space on both sides to set instruments down, and an 8" backsplash to protect against water splashes.
- Three sink bowls are included, which feature a unique, "no drip" sink lip, pitched bottoms for proper drainage, and permanently etched gallon and liter calibration inside each bowl for precise measurement.
- Faucets are all sold separately. Options include: a pre-rinse combination faucet, a wall-mounted faucet, and a wall-mounted air, gas and VAC valve with serrated barbed end hose.

Key Product Features Pre-rinsed combination faucet (sold separately) Wall mount Symmetrical design 8" backsplash 18" counter space on faucet with 12' 25" x 18" x 10" sinks allow for for either left or right both sides swing, and 4" full immersion of instrument work flow wrist blade travs handles (sold separately) 10" deep sinks (internal dimensions) Ergonomic, 37" working height F 4 (33" - 40" range for hydraulic units) Permanently etched gallon and Sinks have a pitched bottom liter calibration inside each bowl for for proper drainage and a precise measurement no-drip sink lip

www.blickman.com

Description: We offer three decontamination stations. One sits on two base cabinets, and the other two sit on stainless steel legs that are either stationary or adjustable. All three types of decontamination stations are set to ergonomic working heights. You have the choice of wall-mounted and deck-mounted faucets. All faucets sold separately.

Decontamination Station External Dimensions


Width: 116" x Height: 37" x Depth: 26 1/2" or 30" *33" - 40" adjustable height range for hydraulic unit

Base Cabinet Decontamination Station

Mounting Types	Wall mounted Faucets 26 1/2" Deep	Deck Mounted Faucets 30" Deep
Model #	DSBC26	DSBC30
Item #	162DSBC001	162DSBC000
Shipping Weight	371 lbs.	391 lbs.

Stationary Decontamination Station

Mounting Types	Wall mounted Faucets 26 ¹ / ₂ " Deep	Deck Mounted Faucets 30" Deep
Model #	DSHB26	DSHB30
Item #	162DSHB001	162DSHB000
Shipping Weight	442 lbs.	462 lbs.

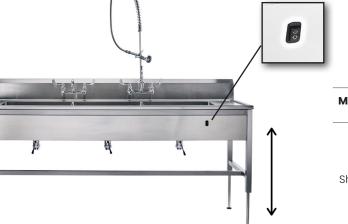

Hick

Hydraulic Decontamination Station

Mounting Types	Wall mounted Faucets 26 ¹ / ₂ " Deeps	Deck Mounted Faucets 30" Deep
Model #	DSHY26	DSHY30
Item #	162DSHY017	162DSHY000
Shipping Weight	442 lbs.	462 lbs.

All decontamination stations include three station sinks.

Sink Dimensions: 25" x 18" x 10"

Options

FAUCET OPTIONS

| All faucets sold separately |

Deck Mounted



Wall Mounted



Name	Pre-rinse Combo Faucet	Swing Faucet	Swivel Gooseneck Faucet	DI/RO Lab Faucet & Lab Hose Fitting Kit
		D	Deck Mount	
Part #	9200601000	9200002214	9243232204	9200586002
Shipping Weight	13 lbs.	5 lbs.	5 lbs.	5 lbs.
		N	Wall Mount	
Part #	9200133010	9200231020	9243232414	9204250KIT
Shipping Weight	13 lbs.	5 lbs.	5 lbs.	2 lbs.

Note: The DI/RO Lab Faucet is not for water, but for gas, air, or vacuum hose purposes

Accessories



Back Drain

Cleaning	
Spray Gun	

Valve

Name	Eyewash	Overflow Back Drain	Cleaning Spray Gun + Tabletop Rosette	Thermo Valve
Part #	9201849LHL	9201074150	9200107006 + 9200800017	9203600LF0
Shipping Weight	8 lbs.	6 lbs.	2 lbs. + 1 lb.	7 lbs.

Note: The Cleaning Spray Gun comes with 1 tip. 8 different sized replacement tips can be ordered: **Part #:** 9200107006

Pegboard

Optional add-on by request



Construction Specifications

	Construction Specifications
Standard Dimensions	External Dimensions: Width: 116" x Depth: 26 ¹ / ₂ or 30" x Height: 37" (33" - 40" hydraulic unit) Internal Sink Bowl Dimensions: Width: 25" x Height: 18" x Depth: 10"
Materials	304 stainless steel
Gauge	Bowl: 14GA Sink Top and Channel: 16GA DSHY Legs: 16GA; Cross Bracing Square Tube: 18GA DSHB Legs and Cross Bracing Square Tube
Welding	All welding by heliarc method and exposed welds are ground smooth
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish
Electrical	Hydraulic Unit: 120V AC
Casters	N/A
Notes	 Ergonomic working height for reduced back stress, symmetrical design for either left or right work flow
	 18" counter space on both sides to set instruments down
	 8" backsplash to protect against water splashes
	• Three sink bowls which feature a unique, "no drip" sink lip, pitched bottoms for proper drainage, and permanently etched gallon and liter calibration inside each bowl for precise measurement
	All faucets sold separately

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Toll Free: (800) 247-5070 Local: (973) 330-0557



www.blickman.com info@blickman.com